


Tuttnauer

LEADING IN
STERILIZATION
SINCE

1925

T-Top Autoclave

Efficiency with maximum reliability

The only autoclave upgradeable to pre/post vacuum cycles. A cost effective & durable closed door drying solution for medical & healthcare facilities.



**Designed to Meet Your Healthcare
Clinics' Everyday Sterilization Needs**

Get it from your authorized distributor: www.labequip.com

Labequip

ask@labequip.com

T-Top 10|11 Autoclave

LEADING IN
STERILIZATION
SINCE

1925

- ◆ Class S sterilizer (Closed door drying)
- ◆ Upgradable to pre/post vacuum (Class B)
- ◆ Clean water with every cycle! Separate water reservoirs for clean /used water
- ◆ T-Top 11 comes with automatic drain and fill kit included
- ◆ **Best price for TOP performance**
- ◆ T-Connect: Advanced remote diagnostics
- ◆ Large chamber capacity



P/N: T-Top 10/11 110V



Innovation

Best in class technology

- ◆ Time saver setting to set tests and sterilization cycles to start at any time
- ◆ Programming and customization functions
- ◆ Specially designed user interface to speed up and optimize daily workload



Plug and Play

No dedicated installation connections required

- ◆ Multicolour user-friendly touch screen with intuitive icons
- ◆ Compact and ergonomic, designed to fit into your space
- ◆ Easy - just start and run programs



Traceability & Connectivity

Advanced Technology that provides real-time data intelligence

- ◆ 2 USB Ports
- ◆ Wi-Fi connection & T-Connect software
- ◆ Remotely access cycle history, get notifications and alerts.
- ◆ Paperless documentation



Reliable Performance

Quality product that will do all you need and expect

- ◆ Easy to operate, maintain and use
- ◆ Internal memory - storing up to 500 cycles
- ◆ Advanced technology to deliver longevity and reliability

Cloud based paperless software

T-Top devices feature the exclusive and technologically advanced T-Connect, a cloud based solution for your sterilization center.




T-Top Sterilization Cycle Programs

Class S	Icon	Cycle program	Load type	T-Top 10" cycle times including drying	T-Top 11" cycle times including drying
		Unwrapped instruments 270°F	Unwrapped Instruments	21 min	37 min
		Wrapped Pouches 270°F	Wrapped Pouches	47 min	59 min
		Unwrapped Delicate 250°F	Unwrapped Delicate instruments	34 min	46 min
		Handpieces 270°F	Handpieces	50 min	70 min
Class B	Icon	Cycle program	Load type	T-Top 10" cycle times including drying	T-Top 11" cycle times including drying
		Unwrapped instruments 273°F	Unwrapped Instruments	24 min	38 min
		Wrapped Pouches 273°F	Wrapped Pouches	52 min	66 min
		Wrapped Pouches 250°F	Wrapped Delicate instruments	73 min	80 min
		Unwrapped Delicate 250°F	Unwrapped Delicate instruments	34 min	46 min

* Cycle times above are an average for a standard load and could vary based on load weight, type, and external conditions

T-Top Technical Specifications

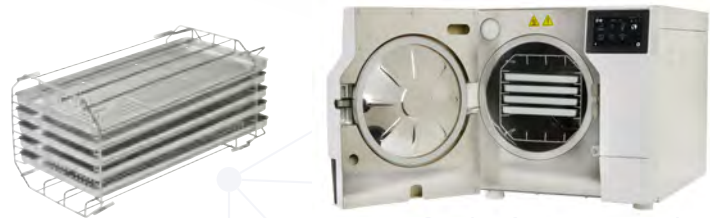

	T-Top 10"	T-Top 11"
Chamber size	10 (W) x 18 (D)	11 (W) x 18 (D)
Volume	21L	27L
Overall dimensions	23" (D) x 18.1" (W) x 18" (H)	23" x 19.5" x 17.5"
Max # of trays	4 Trays	5 Trays
Weight	99 lbs.	115 lbs.
Water filling	Manual	Manual and Automatic

Tuttnauer USA Co.
345 Oser Ave.
Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850
Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com

Visit Our Website:
www.tuttnauer.com/us

